


Potassium Clavulanate

クラブラン酸カリウム

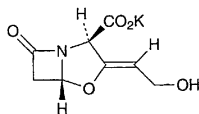

$C_8H_8KNO_5$: 237.25

Monopotassium (2*R*,5*R*)-3-[(1*Z*)-2-hydroxyethylidene]-7-oxo-4-oxa-1-azabicyclo[3.2.0]heptane-2-carboxylate [61177-45-5]

Potassium Clavulanate contains not less than 755 μ g (potency) per mg, calculated on the anhydrous basis. The potency of Potassium Clavulanate is expressed as mass (potency) of clavularic acid ($C_8H_9NO_5$; 199.16).

Description Potassium Clavulanate occurs as a white to light yellowish white, crystalline powder.

It is very soluble in water, soluble in methanol, and slightly soluble in ethanol (95).

It is hygroscopic.

Identification (1) To 1 mL of a solution of Potassium Clavulanate (1 in 50,000) add 5 mL of imidazole TS, warm in a water bath at 30°C for 12 minutes. After cooling, determine the absorption spectrum of this solution as directed under the Ultraviolet-visible Spectrophotometry: it exhibits a maximum between 311 nm and 315 nm.

(2) Determine the infrared absorption spectrum of Potassium Clavulanate as directed in the potassium bromide disk method under the Infrared Spectrophotometry, and compare the spectrum with the Reference Spectrum: both spectra exhibit similar intensities of absorption at the same wave numbers.

(3) Potassium Clavulanate responds to the Qualitative Test (1) for potassium salt.

Optical rotation $[\alpha]_D^{20}$: +55 – +60° (0.5 g calculated on the anhydrous basis, water, 50 mL, 100 mm).

pH Dissolve 0.1 g of Potassium Clavulanate in 10 mL of water: the pH of the solution is between 6.0 and 8.0.

Purity (1) Heavy metals—Proceed with 2.0 g of Potassium Clavulanate according to Method 2, and perform the test. Prepare the control solution with 4.0 mL of Standard Lead Solution (not more than 20 ppm).

(2) Arsenic—Prepare the test solution with 2.0 g of Potassium Clavulanate according to Method 3, and perform the test using Apparatus B (not more than 2 ppm).

Water Not more than 1.5% (5 g, volumetric titration, direct titration).

Assay Weigh accurately an amount of Potassium Clavulanate and Lithium Clavulanate Reference Standard, equivalent to about 12.5 mg (potency), dissolve each in 30 mL of water, add exactly 5 mL of the internal standard solution and water to make 50 mL, and use these solutions as the sample solution and the standard solution, respectively. Perform the test with 5 μ L each of the sample solution and

the standard solution as directed under the Liquid Chromatography according to the following conditions, and calculate the ratios, Q_T and Q_S , of the peak area of clavularic acid to that of the internal standard.

Amount [μ g (potency)] of clavularic acid ($C_8H_9NO_5$)
= amount [mg (potency)] of Lithium Clavulanate

$$\text{Reference Standard} \times \frac{Q_T}{Q_S} \times 1000$$

Internal standard solution—Dissolve 0.3 g of sulfanilamide in 30 mL of methanol, and add water to make 100 mL.

Operating conditions—

Detector: An ultraviolet absorption photometer (wavelength: 230 nm).

Column: A stainless steel column 4.6 mm in inside diameter and 25 cm in length, packed with octadecylsilanized silica gel for liquid chromatography (5 μ m in particle diameter).

Column temperature: A constant temperature of about 25°C.

Mobile phase: Dissolve 1.36 g of sodium acetate trihydrate in 900 mL of water, adjust to pH 4.5 with diluted acetic acid (31) (2 in 5), and add 30 mL of methanol and water to make 1000 mL.

Flow rate: Adjust the flow rate so that the retention time of clavularic acid is about 6 minutes.

System suitability—


System performance: When the procedure is run with 5 μ L of the standard solution under the above operating conditions, clavularic acid and the internal standard are eluted in this order with the resolution between these peaks being not less than 4.

System repeatability: When the test is repeated 6 times with 5 μ L of the standard solution under the above operating conditions, the relative standard deviation of the ratios of the peak area of clavularic acid to that of the internal standard is not more than 1.0%.

Containers and storage Containers—Tight containers.

Potassium Guaiacolsulfonate

ゲアヤコールスルホン酸カリウム

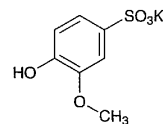

$C_7H_7KO_5S$: 242.29

Monopotassium 4-hydroxy-3-methoxybenzenesulfonate [1321-14-8]

Potassium Guaiacolsulfonate contains not less than 98.5% of $C_7H_7KO_5S$, calculated on the anhydrous basis.

Description Potassium Guaiacolsulfonate occurs as white crystals or crystalline powder. It is odorless or has a slight, characteristic odor and a slightly bitter taste.

It is freely soluble in water and in formic acid, soluble in

methanol, and practically insoluble in ethanol (95), in acetic anhydride and in diethyl ether.

Identification (1) To 10 mL of a solution of Potassium Guaiacolsulfonate (1 in 100) add 2 drops of iron (III) chloride TS: a blue-purple color develops.

(2) Dissolve 0.25 g of Potassium Guaiacolsulfonate in water to make 500 mL, and to 10 mL of this solution add phosphate buffer solution, pH 7.0, to make 100 mL. Determine the absorption spectrum of this solution as directed under the Ultraviolet-visible Spectrophotometry, and compare the spectrum with the Reference Spectrum: both spectra exhibit similar intensities of absorption at the same wavelengths.

(3) A solution of Potassium Guaiacolsulfonate (1 in 10) responds to Qualitative Tests for potassium salt.

pH Dissolve 1.0 g of Potassium Guaiacolsulfonate in 20 mL of water: the pH of the solution is between 4.0 and 5.5.

Purity (1) Clarity and color of solution—Dissolve 1.0 g of Potassium Guaiacolsulfonate in 20 mL of water: the solution is clear and colorless.

(2) Sulfate—Perform the test with 0.8 g of Potassium Guaiacolsulfonate. Prepare the control solution with 0.50 mL of 0.005 mol/L sulfuric acid VS (not more than 0.030%).

(3) Heavy metals—Proceed with 1.0 g of Potassium Guaiacolsulfonate according to Method 1, and perform the test. Prepare the control solution with 2.0 mL of Standard Lead Solution (not more than 20 ppm).

(4) Arsenic—Prepare the test solution with 1.0 g of Potassium Guaiacolsulfonate according to Method 1, and perform the test using Apparatus B (not more than 2 ppm).

(5) Related substances—Dissolve 0.20 g of Potassium Guaiacolsulfonate in 200 mL of mobile phase, and use this solution as the sample solution. Pipet 1 mL of the sample solution, add the mobile phase to make exactly 100 mL, and use this solution as the standard solution. Perform the test with 5 μ L each of the sample solution and the standard solution as directed under the Liquid Chromatography according to the following conditions. Determine each peak area obtained from these solutions by the automatic integration method: the total area of peaks other than the peak of potassium guaiacolsulfonate from the sample solution is not larger than the peak area of potassium guaiacolsulfonate from the standard solution.

Operating conditions—

Detector: An ultraviolet absorption photometer (wavelength: 279 nm).

Column: A stainless steel column 4 mm in inside diameter and 20 to 25 cm in length, packed with dimethylaminopropylsilanized silica gel, (5 to 10 μ m in particle diameter).

Column temperature: A constant temperature of about 30°C.

Mobile phase: A mixture of 0.05 mol/L potassium dihydrogenphosphate VS and methanol (20:1).

Flow rate: Adjust the flow rate so that the retention time of potassium guaiacolsulfonate is about 10 minutes.

Selection of column: Weigh 0.050 g each of potassium guaiacolsulfonate and guaiacol, and dissolve in 50 mL of the mobile phase. Proceed with 5 μ L of this solution under the above operating conditions, and calculate the resolution. Use a column giving elution of guaiacol and potassium

guaiacolsulfonate in this order with the resolution of these peaks being not less than 4.

Detection sensitivity: Adjust the sensitivity so that the peak height of potassium guaiacolsulfonate from 5 μ L of the standard solution is not less than 10 mm.

Time span of measurement: About twice as long as the retention time of potassium guaiacolsulfonate.

Water 3.0 – 4.5% (0.3 g, direct titration).

Assay Weigh accurately about 0.3 g of Potassium Guaiacolsulfonate, dissolve in 2.0 mL of formic acid, add 50 mL of acetic anhydride, and titrate with 0.1 mol/L perchloric acid VS (potentiometric titration). Perform a blank determination, and make any necessary correction.

Each mL of 0.1 mol/L perchloric acid VS
= 24.229 mg of C₇H₇KO₅S

Containers and storage Containers—Well-closed containers.

Storage—Light-resistant.

Potassium Iodide

ヨウ化カリウム

KI: 166.00

Potassium Iodide, when dried, contains not less than 99.0% of KI.

Description Potassium Iodide occurs as colorless or white crystals, or a white crystalline powder.

It is very soluble in water, soluble in ethanol (95), and practically insoluble in diethyl ether.

It is slightly deliquescent in moist air.

Identification A solution of Potassium Iodide (1 in 20) responds to the Qualitative Tests for potassium salt and for iodide.

Purity (1) Clarity and color of solution—Dissolve 1.0 g of Potassium Iodide in 2 mL of water: the solution is clear and colorless.

(2) Alkali—Dissolve 1.0 g of Potassium Iodide in 10 mL of freshly boiled and cooled water, and add 0.50 mL of 0.005 mol/L sulfuric acid and 1 drop of phenolphthalein TS: no color develops.

(3) Chloride, bromide and thiosulfate—Dissolve 0.20 g of Potassium Iodide in 5 mL of ammonia TS, add 15.0 mL of 0.1 mol/L silver nitrate VS, shake for 2 to 3 minutes, and filter. To 10 mL of the filtrate, add 15 mL of dilute nitric acid: no brown color develops. The solution has no more turbidity than that of the following control solution.

Control solution: To 0.30 mL of 0.01 mol/L hydrochloric acid VS add 2.5 mL of ammonia TS, and 7.5 mL of 0.1 mol/L silver nitrate VS and 15 mL of dilute nitric acid.

(4) Nitrate, nitrite and ammonium—Place 1.0 g of Potassium Iodide in a 40-mL test tube, and add 5 mL of water, 5 mL of sodium hydroxide TS and 0.2 g of aluminum wire. Insert the absorbent cotton in the mouth of the test tube, and place a piece of moistened red litmus paper on it. Heat the test tube carefully on a water bath for 15 minutes: